


Quick Guide











Gmate® Origin



Side Buttons(Up/Down Button)
are used to recall the stored results
and adjust the parameters in setting
mode.

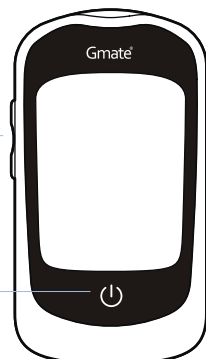

Power Button
is used to meter on/off and
change the meter modes.

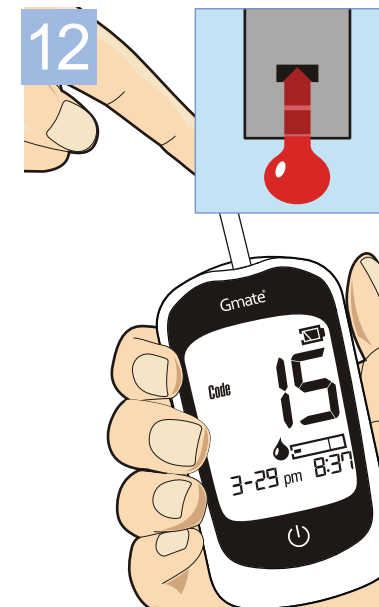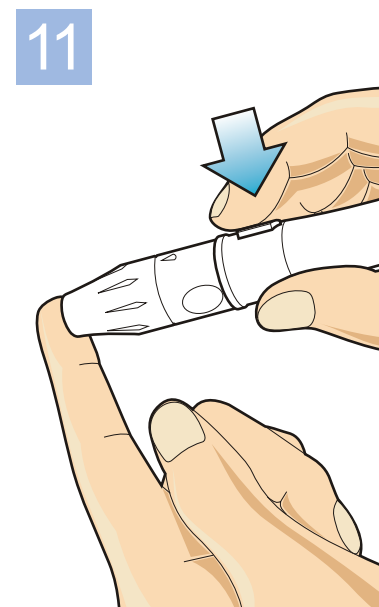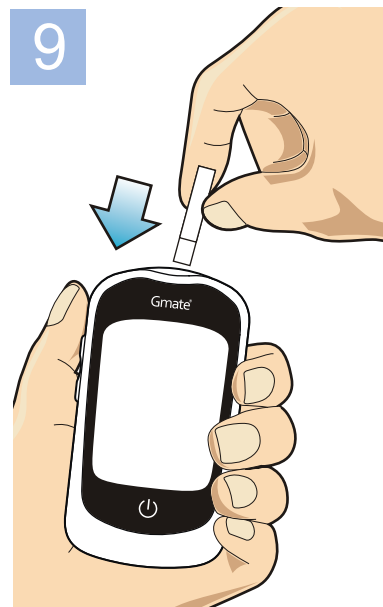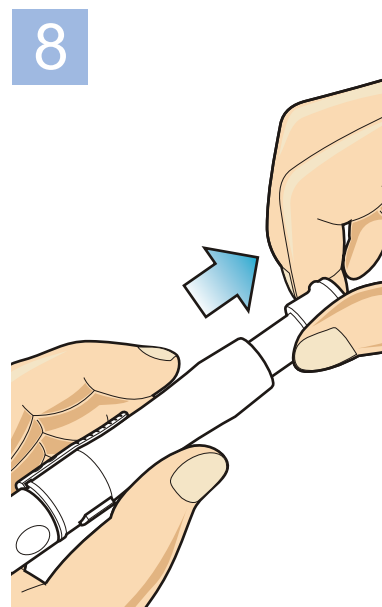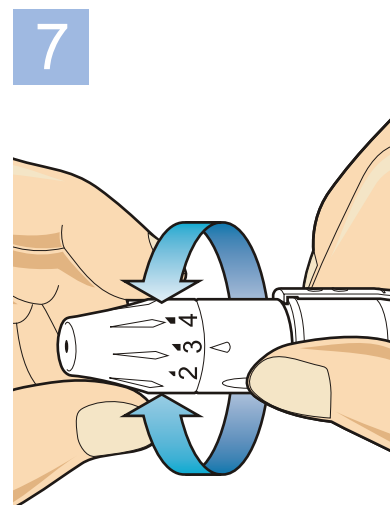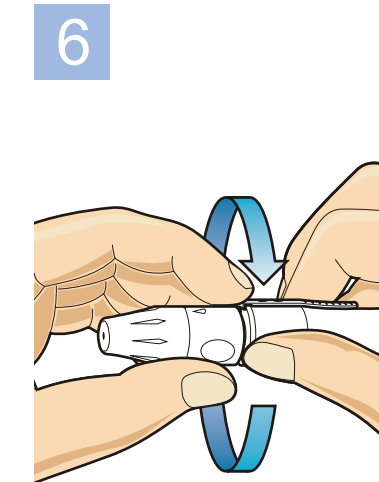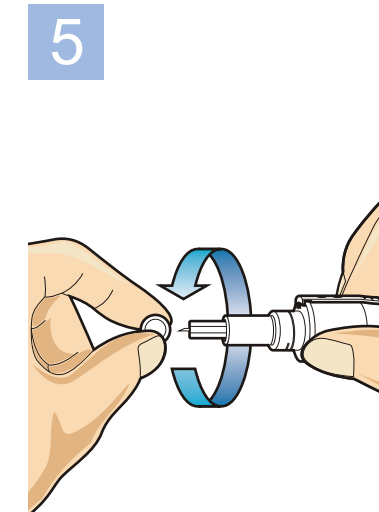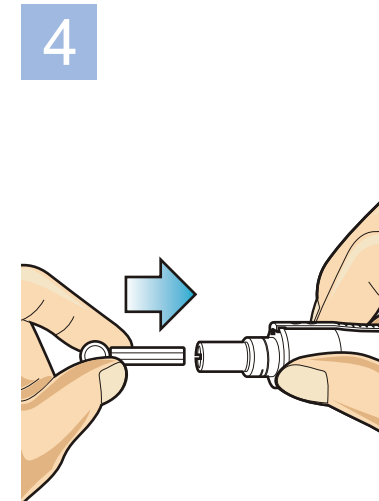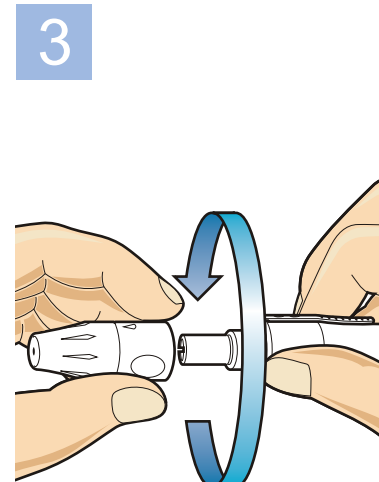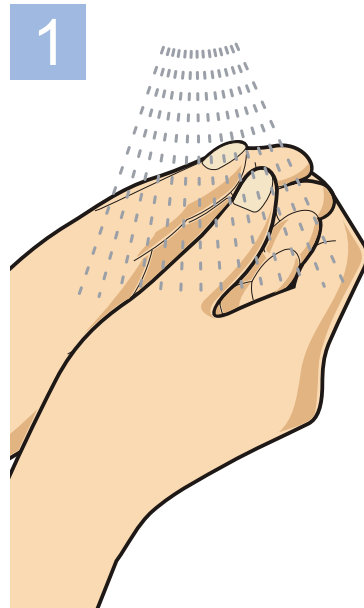

- E01** When inserting a used test strip.
- E02** When the test strip may be inserted incorrectly or have a problem.
- E03** When the meter has problem.
- E04** When environmental temperature is BELOW or ABOVE system operation range, (10~32°C)
- E05** Error message that may appear during your blood glucose testing. Please check "Testing Your Blood Glucose" page in the User's Guide and repeat the test with a new strip.